

# Sapphire SF NX

# CW Free Space and Fiber Pigtailed Single-Frequency Lasers

Sapphire SF NX is a series of compact CW visible single-frequency lasers based on Coherent's unique OPSL (Optically Pumped Semiconductor Laser) technology. Sapphire SF NX lasers come with an ultranarrow linewidth of <1.5 MHz, a high quality diffraction-limited beam with excellent pointing stability, high power stability and low noise.

Sapphire lasers are manufactured in cleanrooms using Coherent's patented PermAlign™ technology for optimal aligning and solder-bonding the optics. Sapphire SF NX lasers come with a flexible interface concept: Analog, RS-232 or USB – it is up the user to select the appropriate communication channel.

Sapphire SF NX lasers are intended for applications that need narrow and ultra-narrow linewidth light such as Raman spectroscopy, interferometry, holography, metrology, and inspection.



## **FEATURES & BENEFITS**

- Ultra-narrow linewidth
- Outstanding power and wavelength stability
- Output power
  - Up to 150 mW at 488 nm / 100 mW fiber pigtailed
  - Up to 200 mW at 532 nm / 150 mW fiber pigtailed
- · Ultra-low noise
- Superior beam quality
- PermAlign technology
- Permanent optimal alignment
- Unsurpassed robust and stable
- · OEM and end-user versions
- Proven Sapphire reliability

#### **APPLICATIONS**

- Raman Spectroscopy
- Interferometry
- Holography
- Metrology
- Inspection



SPECIFICATIONS	Sapphire 488 SF NX	Sapphire 532 SF NX
Wavelength <sup>1</sup> (nm)	488.0	532.0
Wavelength Accuracy <sup>2</sup> (nm)	±0.1	
Single-longitudinal Mode, Linewidth (MHz)	<1.5	
Output Power <sup>3</sup> (mW)	20, 50, 75, 100, 150	20, 50, 75, 100, 150, 200
Spatial Mode	TEM <sub>00</sub> , 1	$M^2 < 1.1$
Beam Asymmetry	0.9 to 1.1	
Beam Diameter at 1/e <sup>2</sup> (mm)	0.70 ±0.05	
Beam Divergence (mrad)	<1.3	
Pointing Stability (µrad) (over 2 hours after warm-up and ±3°C)	<30	
Noise (%)		
20 Hz to 2 MHz, rms	<0.25	
20 Hz to 20 kHz, peak-to-peak	<1	
Long-term Power Stability (%) (2 hours, ±3°C)	<2	
Warm-up Time (minutes)	<5	
Polarization Ratio	>100:1, vertical	
Static Alignment Tolerances <sup>4</sup>		
Beam Position (mm)	±0.25	
Beam Angle (mrad)	±2.5	
Beam Waist Position with respect to Exit Window	±20	005
UTILITY AND ENVIRONMENTAL REQUIR	EMENTS	
Interfacing	Analog, RS-232, USB	
Operating Voltage <sup>6</sup> (VDC)		
Maximum Rated	+10.8 to 15.0	
Nominal	+12.0 to 13.2	
Power Consumption (W)	<60	
Max. Laser Head Baseplate Temperature <sup>7</sup>	+50°C (122°F)	
Max. Heat Dissipation of Head (W)	25 (baseplate at 50°C)	
Ambient Temperature		
Operating Condition	10 to 40°C (50 to 104°F) non-condensing	
Non-Operating Condition	-30 to 60°C (-22 to 140°F)	
Shock Tolerance (6 ms)	7 g laterally, 15 g vertically	
Dimensions (L x W x H)	125 70 24	(4.0
Laser Head Controller	125 x 70 x 34 mm (4.9 x 2.8 x 1.3 in.)	
Heat Sink (optional)	123.3 x 76.2 x 38.1 mm (4.9 x 3.0 x 1.5 in.) 200 x 80 x 50 mm (7.9 x 3.2 x 2 in.)	
DC Power Supply (optional)	171 x 104 x 55 mm (6.7 x 4.1 x 2.2 in.)	
Cable — Laser Head to Controller	1.52 m (5.0 ft.)	
Weights	1.52	,
Laser Head	0.35 kg (0.77 lbs.)	
Controller	0.25 kg (0.55 lbs.)	
Heat Sink (optional)	0.75 kg (1.65 lbs.)	
DC Power Supply (optional)	0.95 kg (incl. line cable) (2.1 lbs.)	
Packaged System (head+controller+cable+manual)	1.7 kg (3.7 lbs.)	
Cable — Laser Head to Controller	0.3 kg (0.66 lbs.)	

<sup>1</sup> Measured in air.




Laser-to-laser tolerance and wavelength tolerance over guaranteed lifetime.

Laser-to-baser torerainte aniu wavelengtin Toleraince over guaranteed lifetime.
 Output power is adjustable via analog or digital interface from 10% to 110%. Specifications are valid for 100% power.
 Static alignment tolerances are relative to the right bottom edge (in beam direction).
 200 mm is ~25% of Rayleigh Range at 488 mm and ~30% of Rayleigh Range at 532 mm.
 If user-supplied, the DC power supply has to meet the following requirements: Power >60W; ripple <5% peak-to-peak; line regulation <0.5%.</li>
 With factory-provided or other adequate heat sink.

### **MECHANICAL SPECIFICATIONS**

# Sapphire SF NX

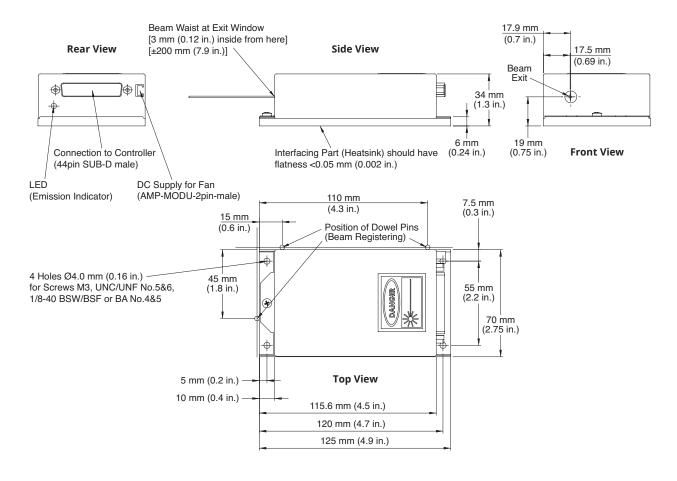




SPECIFICATIONS	Sapphire 488 SF FP NX	Sapphire 532 SF FP NX
Wavelength¹ (nm)	488.0	532.0
Wavelength Accuracy <sup>2</sup> (nm)	±0.1	
Single-longitudinal Mode, Linewidth (MHz)	<1.5	
Output Power at Fiber Exit <sup>3</sup> (mW)	20, 50, 100	20, 50, 100, 150
Fiber Type	SM-PM <sup>4</sup>	
Fiber Length (m)	1	
Fiber Numerical Aperture (NA) (1/e <sup>2</sup> )	0.06	
Fiber Output	FC/APC; 8° angled <sup>5</sup>	
Spatial Mode	TEM <sub>00</sub> , M <sup>2</sup> <1.1	
Beam Asymmetry	<1:1.1	
Noise (%) 20 Hz to 2 MHz, rms 20 Hz to 20 kHz, peak-to-peak	<0.25 <1	
Long-term Power Stability (%) (2 hours, ±3°C)	<2	
Warm-up Time (minutes)	<5	
Polarization Ratio	>100:1, line	ar, vertical
UTILITY AND ENVIRONMENTAL REQUIREME		
Interfacing	Analog, RS-232, USB	
Operating Voltage <sup>6</sup> (VDC)  Maximum Rated  Nominal	+10.8 to 15.0 +12.0 to 13.2	
Power Consumption (W)	<60	
Max. Laser Head Baseplate Temperature <sup>7</sup>	+50°C (122°F)	
Max. Heat Dissipation of Head (W)	25 (baseplate at 50°C)	
Ambient Temperature Operating Condition Non-Operating Condition	10 to 40°C (50 to 104°F) non-condensing -30 to 60°C (-22 to 140°F) 7 g laterally, 15 g vertically	
Shock Tolerance (6 ms)	/ g laterally, 15	g vertically
Dimensions (L x W x H)  Laser Head <sup>8</sup> Controller  Heat Sink (optional)  DC Power Supply (optional)  Cable — Laser Head to Controller	125 x 70 x 34 mm (4.9 x 2.8 x 1.3 in.) 123.3 x 76.2 x 38.1 mm (4.9 x 3.0 x 1.5 in.) 200 x 80 x 50 mm (7.9 x 3.2 x 2 in.) 171 x 104 x 55 mm (6.7 x 4.1 x 2.2 in.) 1.52 m (5.0 ft.)	
Weights Laser Head <sup>8</sup> Controller Heat Sink (optional) DC Power Supply (optional) Packaged System (head+controller+cable+manual) Cable — Laser Head to Controller	0.35 kg (0.77 lbs.) 0.25 kg (0.55 lbs.) 0.75 kg (1.65 lbs.) 0.95 kg (incl. line cable) (2.1 lbs.) 1.7 kg (3.7 lbs.) 0.3 kg (0.66 lbs.)	

- Measured in air.
   Laser-to-laser tolerance and wavelength tolerance over guaranteed lifetime.
   Output power is adjustable via analog or digital interface from 10% to 110%. Specifications are valid for 100% power.
   Single-mode, polarization maintaining fiber, bending radius min. 50 mm.
   Fiber FCAPC connector output not compatible for patchcord-to-patchcord connection.
   If user-supplied, the DC power supply has to meet the following requirements: power >60W; ripple <5% peak-to-peak; line regulation <0.5%.</li>
   With factory-provided or other adequate heat sink.
   Dimposions respectively weighth without fiber virtal narr
- 8 Dimensions respectively weight without fiber pigtail part.




# **MECHANICAL SPECIFICATIONS**

# Sapphire SF FP NX

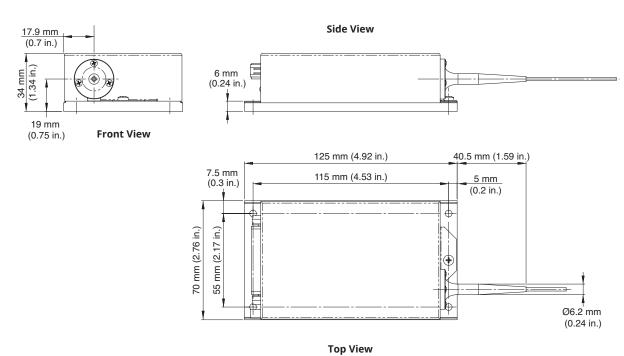

#### MECHANICAL SPECIFICATIONS

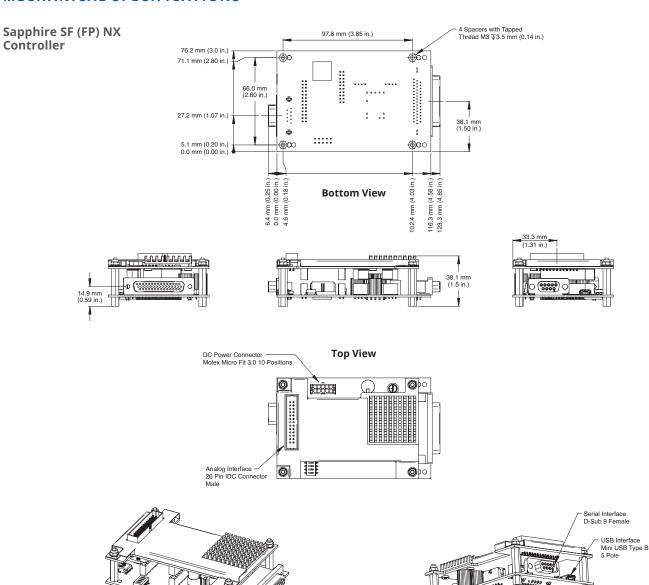

Laser Head D-Sub 44 Female



Coherent, Inc.,

5100 Patrick Henry Drive Santa Clara, CA 95054

p. (800) 527-3786 | (408) 764-4983

f. (408) 764-4646

### tech.sales@coherent.com www.coherent.com

VISIBLE AND INVISIBLE LASER RADIATION.

DIRECT OR SCATTERED RADIATION.

LASE IN LASER PRODUCT.

AMEL. Power

1. 51 to 0.50 pm 400 www.

0.90 to 1.00 pm 50 mW

**( E** 



Coherent follows a policy of continuous product improvement. Specifications are subject to change without notice. Coherent's scientific and industrial lasers are certified to comply with the Federal Regulations (21 CFR Subchapter J) as administered by the Center for Devices and Radiological Health on all systems ordered for shipment after August 2, 1976.